

Elbow type 1
For common tube (single limb)



Elbow type 2
For professional medical breathing circuit (dual limb)

Soft and Comfortable



Newly upgraded single layer liquid silicone process; Anti-suffocation safety valve; safer to wear and more comfortable; 360° rotatable detachable elbow

Easy to use



Headgear buckle can be quickly removed for easy wearing and removal; Disassembly and cleaning are easier.

Safe to use



Medical grade odorless silicagel, safe material, feel free to use; Soft and skin-friendly headgear, passed BIOC test, suitable for sensitive skin.

Wearing Instructions

Step 1



Put the headgear through the mask from both sides.

Step 2



Rightly wear the mask on head and fix it

Step 3




Tie up the lower parts of the headgear to other two buckles

Step 4



Adjust the position of headgear and forehead pad to make the mask perfectly fit.

Dimension and Accessory Replacement

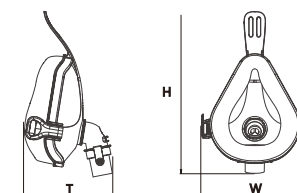

Dimension Size	H(mm)	W(mm)	T(mm)
Large	186.8	123	99
Medium	181.2	122	99.1
Small	163.4	114	96.8

Accessory Replacement Suggestion	Silicon cushion >> Every 6 months	Headgear >> Every 6 months	Mask frame >> Every 12-18 months
----------------------------------	-----------------------------------	----------------------------	----------------------------------

EaseFit FM3-NV Face Mask

S M L

Technology Serves Human Health

24 hours service hotline // +86-731-82564299 / +86-19807315856 (WeChat & WhatsApp)

Email : sales@beyondmedtech.com URL: www.beyondmedtech.com

Add:Beyond Zone, Zhixian Rd, Xueshi Street, Yuelu District, 410208, Changsha, China

Version number: CY-02-57-02 A/0-002



Beyond Website