




EaseFit FMI-NV Face Mask

S M L

Technology Serves Human Health



Product Feature

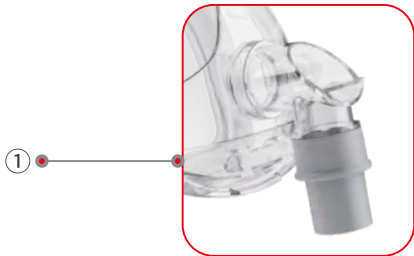

Elbow type 1
For common tube (single limb)



Elbow type 2
For professional medical breathing circuit (dual limb)

Adjustable Height



2 levels of height are available.

Forehead Stepless Knob



Independent invention patent; Stepless knob design; Perfectly suit individual facial contour.

Double Cushions



Medical liquid silicone; A comfortable seal is achieved by dual-wall design; Easy to clean and reassemble.

Arc-buckle



Curve buckle structure enables flexible rotation; Easy to attach, easy to adjust.

Wearing Instructions

Step 1



Put the headgear through the mask from both sides.

Step 2



Put the mask on the head in good position.

Step 3



Adjust the headgear and forehead pad to ensure fixation and comfortability.

Dimension and Accessory Replacement



EaseFit FMI-NV (Elbow type 1) EaseFit FMI-NV (Elbow type 2)

Dimension Size	L(mm)		W(mm)	T(mm)
	Elbow type 1	Elbow type 2		
Large	194.9	180.4	122	114.5
Medium	183.9	169.4	118	114.5
Small	175.9	156.9	110	112.5

Accessory Replacement Suggestion

Silicone cushion >> Every 6 months Headgear >> Every 6 months Mask frame >> Every 12-18 months

