

**BYOND**

High Flow Nasal Cannula Oxygen Therapy

# NC Series Nasal Cannula

|| Designed for patients requiring high flow nasal cannula (HFNC) therapy ||



Effective

Comfortable

Easy to use

Safe

Technology Serves Human Health

# NC1 Series Nasal Cannula

## NC1-S

Air outlet spacing :  $12 \pm 3$  mm  
Applicable flow range: 10-60 (L/min)

## NC1-M


Air outlet spacing :  $16 \pm 3$  mm  
Applicable flow range: 10-80 (L/min)

## NC1-L

Air outlet spacing :  $20.4 \pm 3$  mm  
Applicable flow range: 10-80 (L/min)


### Innovative design, accurate natural fit

The oxygen suction mouth of nasal plug adopts elliptical design, the end is cohesive, and fits the natural curve of human nasal cavity. It can effectively avoid nasal damage, ensure smooth gas flow, and provide a more fitting and safer oxygen inhalation experience.

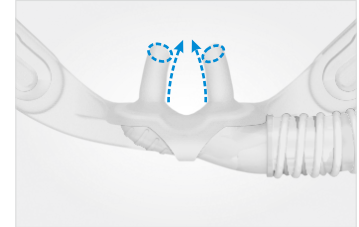


### Ergonomic design, comfortable and light weight to wear

Based on the optimized design of the face model, It reduces the contact area, lowers the friction and pressure between the nose plug and the nasal septum. Even wearing for long periods of time, it significantly improves comfort and provides users a burden-free experience.

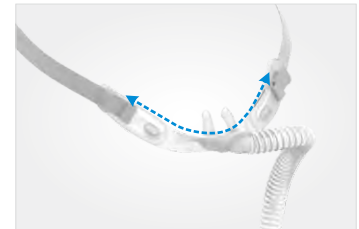

### Stable support, long-lasting and stable gas flow

The elastic headband and high-quality silicone form a stable support structure, which is firm and reliable, reducing the risk of falling off. Combined with the anti-twist and anti-pressure threaded oxygen hose, it effectively avoids airway obstruction and ensures continuous and smooth gas supply.

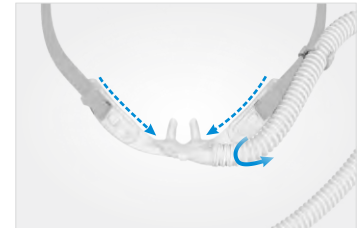

### Medical silicone, skin-friendly and safe protection

Made of medical-grade silicone material, soft and skin-friendly, no smell, through strict biological tests, suitable for people with sensitive skin. Providing users with safe and comfortable use guarantees.



### Compatible with dual-tube and could use with a variety of devices

Supports two connector sizes of 15mm and 22mm inner diameter, compatible with a variety of breathing tubes. Flexible design adapts to different medical equipment and meets the needs of a variety of scenarios, making use more convenient and efficient.



**24<sup>hour</sup> service hotline** // +86-731-82564299 / +86-19807315856 (WeChat & WhatsApp)

Email : [sales@beyondmedtech.com](mailto:sales@beyondmedtech.com) URL: [www.beyondmedtech.com](http://www.beyondmedtech.com)

Add: Beyond Zone, Zhixian Rd, Xueshi Street, Yuelu District, 410208, Changsha, China

Version Number: CY-02-50-02 A/0-002



Beyond Website