

Parameters

Treatment	OSA/COPD		
Model (3.5 Inch)	B18	B17	B19
Work Mode	CPAP, S/T	CPAP,S,T S/T	CPAP,S,T S/T
Screen	black and white screen	3.5-inch TFT color screen	3.5-inch TFT color screen
Pressure Range	4-25cmH ₂ O	4-20cmH ₂ O	4-25cmH ₂ O
SPO ₂ Kit	Optional		
Heating Tube	Optional		
WIFI	Optional		
Data Management Cloud Platform	•		
Ramp	0-60 minutes		
Auto ON/OFF	•		
Auto Leak Compensation	•		
Auto Altitude Adjustment	•		
Noise Level	≤30dBA		
Humidifier	0-5 Level		
Size	280mm*200mm*112mm		
Weight	1KG		

BYOND



Dreamy BIPAP

Technology Serves Human Health

The Dreamy series of BiPAP devices provide customized solutions for adult patients suffering from snoring, sleep apnea, hypopnea syndrome, chronic obstructive pulmonary disease (COPD), respiratory insufficiency, mild to moderate respiratory failure, etc.




Optional Screen




Auto On/Off



Heating And Humidification

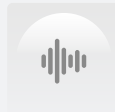

Quiet And Comfortable

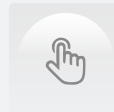

Convenient Operation



◉ Separation Of Main Unit And Humidifier

The separate design of the main unit and humidifier makes it easy to carry, allowing for comfortable respiratory therapy anytime, anywhere.




◉ Humidification Function

Provides 0-5 levels of humidification to adapt to different environments and create a more comfortable breathing experience.




◉ Air leakage Compensation ◉ Ramp

Responding to common leakage issues, it automatically compensates for airflow to prevent pressure reduction and ensure effective therapy.

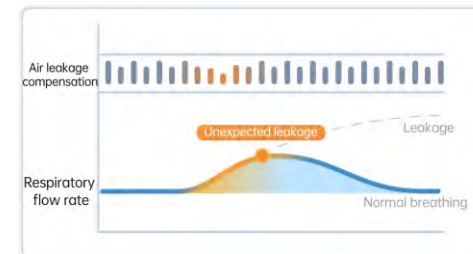


Provides a 0-60 minute ramp function to enhance user tolerance and facilitate falling asleep.

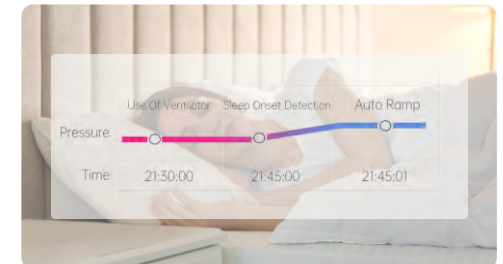

◉ Respirato Synchronization Technique

Automatically follow the user's breathing, accurately judge the trigger and withdrawal time, and quickly adjust to provide the corresponding expiratory pressure and inspiratory pressure, so that the user can breathe more.

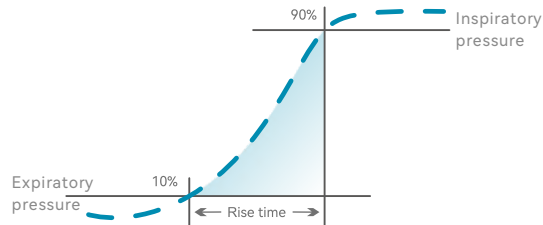

◉ Belex

Triggered by the user's expiratory, the Belex function helps expel bodily waste gases, making breathing easier.

