

Parameters

Treatment	OSA			
Model	C01	C02	AC08	AC09
Work Mode	CPAP	CPAP	CPAP Auto CPAP	CPAP Auto CPAP
Screen	3.5-inch TFT color screen	black and white screen	black and white screen	3.5-inch TFT color screen
Pressure Range	4-20cmH ₂ O	4-20cmH ₂ O	4-20cmH ₂ O	4-20cmH ₂ O
SPO ₂ Kit	Optional			
Heating Tube	Optional			
WIFI	Optional			
Data Management Cloud Platform	•			
Ramp	0-60 minutes			
Auto ON/OFF	•			
Auto Leak Compensation	•			
Auto Altitude Adjustment	•			
Noise Level	≤30dBA			
Humidifier	0-5 Level			
Size	280mm*200mm*112mm			
Weight	1KG			

24^{hour} service hotline // +86-731-82564299 / +86-19807315856 (WeChat & WhatsApp)

Email : sales@beyondmedtech.com URL: www.beyondmedtech.com

Add: Beyond Zone, Zhixian Rd, Xueshi Street, Yuelu District, 410208, Changsha, China

Version number: CY-02-03-01 A/0-002



Beyond Website

BYOND



Dreamy CPAP
Technology Serves Human Health



Compact And Portable

The main unit, including the humidifier, weighs approximately 1KG. The detachable design is convenient for travel or outdoor carrying.



Quiet And Comfortable

≤30dB, low-noise design provides users with a quieter sleep environment.

Detachable Design

Detachable humidifier makes daily use more convenient.




Easy Operation

HD 3.5-inch color screen, clear and simple work page, simple shuttle button operation, one-touch start/stop.



Air leakage Compensation

Responding to common leakage issues, it automatically compensates for airflow to prevent pressure reduction and ensure effective therapy.

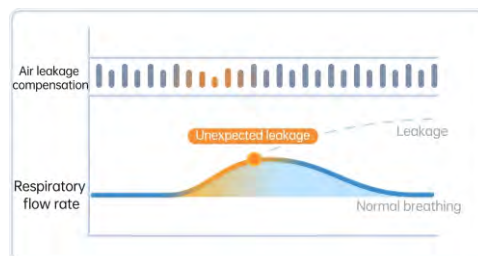

Ramp

The ramp function helps the user fall asleep easily and feel more comfortable.

Data Monitoring

Multiple data monitoring, including pressure and flow waveforms, with 360 days of data storage and 24-hour real-time display.



Flame-Retardant Shell

The new V1-level flame-retardant shell prevents open flames and spontaneous combustion, ensuring safer use.



Intelligent Pressure Regulation

Core technology of pressure and flow triggering, human-machine synchronization, intelligent pressure adjustment, and rapid response.

