

Parameters

Treatment	OSA/COPD		
Model (3.5 Inch)	B-25T	B-30T	B-30P
Work Mode	CPAP,S,T S/T	CPAP,S,T S/T	CPAP,S,T S/T APCV
Pressure Range	4-25cmH ₂ O	4-30cmH ₂ O	4-30cmH ₂ O
SPO ₂ Kit	Optional		
Tidal Volume	—	—	200-2000ml
WIFI	Optional		
Data Management Cloud Platform	Optional Wi-Fi enables cloud platform usage		
Ramp	0-60 minutes		
Auto ON/OFF	•		
Auto Leak Compensation	•		
Auto Altitude Adjustment	•		
Noise Level	≤30dBA		
Humidifier	0-5 Level		
Size	280mm*200mm*146mm		
Weight	2.2KG		

24^{hour} service hotline // +86-731-82564299 / +86-19807315856 (WeChat & WhatsApp)

Email : sales@beyondmedtech.com URL: www.beyondmedtech.com

Add: Beyond Zone, Zhixian Rd, Xueshi Street, Yuelu District, 410208, Changsha, China

Version number: CY-02-08-02 A/0-002



Beyond Website


BYOND



ResPlus BIPAP

Technology Serves Human Health

The ResPlus series of BiPAP devices provide customized solutions for adult patients suffering from snoring, sleep apnea, hypopnea syndrome, chronic obstructive pulmonary disease (COPD), respiratory insufficiency, mild to moderate respiratory failure, etc.

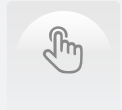

Simple Operation




Heated Humidification



Quiet And Comfortable Use



Tidal Volume Control

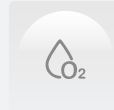


Blood Oxygen Module



Leak Compensation

Real-time monitoring of air leakage provides precise pressure compensation mechanisms to ensure treatment quality and enhance user comfort.

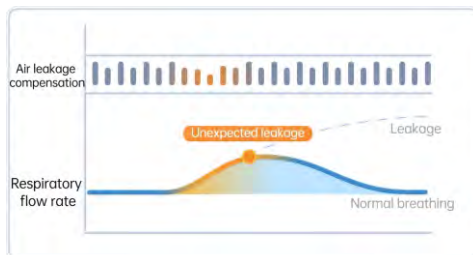



Multiple Modes

A variety of ventilation modes including CPAP, S, AutoB, T, S/T, APCV, etc., to accommodate a wider range of respiratory insufficiencies.




Blood Oxygen Module

Optional blood oxygen module allows users to monitor real-time blood oxygen data, helping them understand the effectiveness of the treatment.

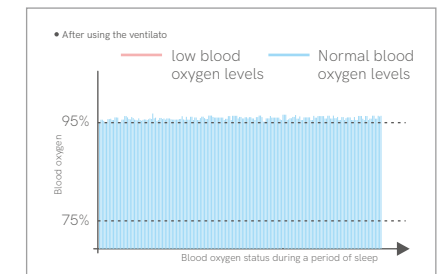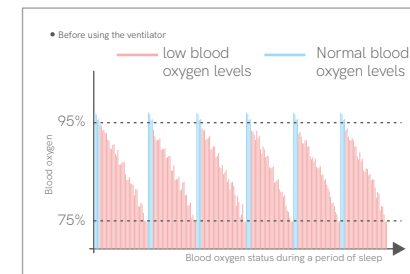

Respirato Synchronization Technique

Automatically follow the user's breathing, accurately judge the trigger and withdrawal time, and quickly adjust to provide the corresponding expiratory pressure and inspiratory pressure, so that the user can breathe more.

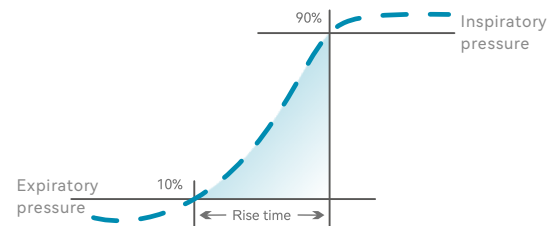

VAF Tidal Volume Control

VAF precise target tidal volume algorithm intelligently detects the required target tidal volume and adjusts it in real-time, ensuring the safety and well-being of the user.

