

Parameters

Treatment	OSA		
Model (3.5 Inch)	C-20C	C-20A	B-25A
Work Mode	CPAP	CPAP APAP	CPAP,S AutoB
Pressure Range	4-20cmH ₂ O	4-20cmH ₂ O	4-25cmH ₂ O
SPO ₂ Kit	Optional		
WIFI	Optional		
Data Management Cloud Platform	Optional Wi-Fi enables cloud platform usage.		
Ramp	0-60 minutes		
Auto ON/OFF	•		
Auto Leak Compensation	•		
Auto Altitude Adjustment	•		
Noise Level	≤30dBA		
Humidifier	0-5 Level		
Size	280mm*200mm*146mm		
Weight	2.2KG		

24^{hour} service hotline // +86-731-82564299 / +86-19807315856 (WeChat & WhatsApp)

Email : sales@beyondmedtech.com URL: www.beyondmedtech.com

Add: Beyond Zone, Zhixian Rd, Xueshi Street, Yuelu District, 410208, Changsha, China

Version number: CY-02-08-01 A/0-002



Beyond Website

BYOND



ResPlus CPAP
Technology Serves Human Health



Compact And Portable Design

The host weight 2.2 kgs with the humidifier and 1.5 kgs without it; The detachable design makes it convenient for travel or portability.



Quiet And Comfortable Operation

≤30 dB, low-noise design, providing users with a quieter sleep environment.



Stable Performance

High-quality sensors ensure stable operation, durability, and fast data responsiveness.

Detachable Design

Separate design of the host and humidifier ensures easy portability, allowing comfortable respiratory therapy anytime, anywhere.



User-Friendly Operation

Equipped with a high-definition 3.5-inch color display with clear and intuitive interface. The dial button enables simple operation with one-touch start/stop function.



Blood Oxygen Module

Optional blood oxygen module allows users to monitor real-time blood oxygen data, helping them understand the effectiveness of the treatment.




Data Monitoring

Comprehensive monitoring data including pressure and flow waveforms, with 360 days of data storage and real-time 24-hour display.



Inspiratory Triggering

Highly sensitive inspiratory triggering technology accurately responds to each breath, reducing breathing resistance.

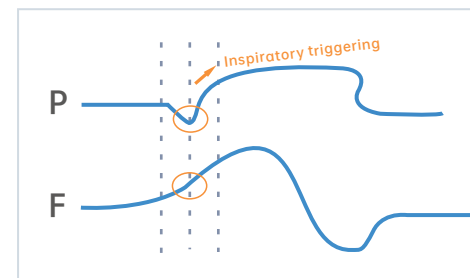

Leak Compensation

Automatic response to common air leakage issues with flow compensation to maintain pressure and ensure treatment effectiveness.

