

Parameters

Model	OG-3A	OG-3C	OG-5A	OG-5C
Oxygen Flow Rate	3	3	5	5
Screen Size	7 Inch LED	7 Inch LED	7 Inch LED	7 Inch LED
Real-Time Oxygen Concentration Display	●	●	●	●
Voice Broadcast Function	●	●	●	●
Nebulization Function	—	●	—	●
Remote Control Function	—	●	—	●
Timing Function	●	●	●	●
Universal Wheels	●	●	●	●
Touch Screen	●	●	●	●
Scheduled Timing Function	30min-24h			
Weight	18±3		19±3	
Noise Level	≤42			



Beyond Medical Molecular Sieve Oxygen Generator



Oxygen inhalation helps maintain good health.

It is the first choice for individuals with moderate to severe hypoxia. High-flow oxygen supply improves blood oxygen saturation. The 5L oxygen concentrator is suitable for oxygen inhalation across various symptoms and age groups.

The Beyond Medical Molecular Sieve Oxygen Generator OG models (3L/5L/min) fully comply with national medical oxygen therapy standards. They are equipped with high-powered oil-free compressors and a robust core, ensuring 7*24h uninterrupted oxygen supply!

* Select the best from the best of the best

Using high-quality imported French lithium sieve (Arkema CECA) and adopting the highly efficient SXSDM type with superior dynamic strength as the OG molecular sieve material solution, achieving rapid PSA oxygen-nitrogen adsorption and separation. (The SXSDM type has been widely used in major oxygen concentrator manufacturers in Europe and America with an outstanding reputation!)

CECA Lithium Sieve· SXSDM Type Solution

High consistency in particle size, smooth and well-rounded, excellent uniformity, and consistent color

Common Lithium Sieve Solutions on the Market

Significant particle size deviations, square mouth cracks, poor uniformity, and inconsistent color

Utilizing high-frequency vibration filling technology to reduce ineffective cavities in the tank, customized intelligent equipment and production lines, along with a strict quality control process for filling.

- ◆ Anti-backflow external structural design helps reduce the risk of molecular sieve moisture absorption.
- ◆ Anti-backflow internal one-way valve design prevents water vapor backflow during operation.
- ◆ Anti-backflow internal normally closed solenoid valve design prevents water vapor backflow when not powered.



* ≤42dB low noise oxygen supply

Seven-layer noise reduction technology: air intake silencer, air outlet silencer, low-noise cooling fan, shell soft package noise reduction, compressor shock absorption silencer system, compressor wrapping noise reduction, S air duct design.



* Convenient and quick oxygen therapy

Three scene modes can be switched with one button, and health is no longer based on feelings. Different users can adjust the oxygen flow rate according to the doctor's prescription, and switch between different scene states instantly to meet personalized oxygen production needs.

Comfort mode

Suitable for sitting/reading

Stretch mode


Suitable for yoga/after exercise

Sleep mode

Suitable for rest at night

* Cloud services, Data management

The "Cloud Platform" smart sharing platform is a good assistant for dealers and a good steward for patients. The all-round ecosystem layout enables Beyond to provide users with more complete, smarter and more diverse products and service solutions.

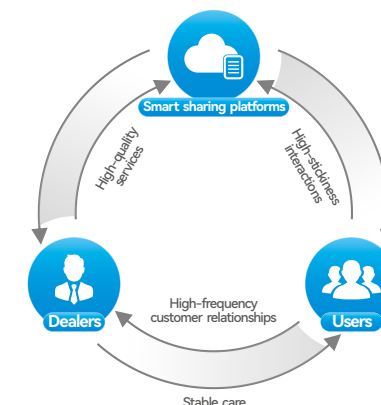

Dealer PC terminal

Through private domain, closed source, exclusive, safe and stable access to Beyond's latest information, customer opportunities, and store operation reports.

User mobile terminal

It can remotely start and stop oxygen therapy instructions (with timing and flow adjustment), monitor the patient's treatment throughout the process, regularly evaluate the treatment effect, remotely send treatment data, and monitor blood oxygen data (later versions will be available), so that problems can be discovered promptly manner and effectively intervened, which can also make the children of older patients feel at ease.

